

KIT Gene Probe Detection Kit

(CW-342)

Intended use

This kit uses orange fluorescein to label KIT orange probe and use green fluorescein to label CEP4 green probe. KIT probe can be combined with the target detection site by in situ hybridization.

Product composition

The kit consists of KIT/CEP4 dual color probe (100µL/Tube).

Storage condition

Keep sealed away from light at $-20^{\circ}\text{C}\pm 5^{\circ}\text{C}$. The product is valid for 12 months. Avoid unnecessary repeated freezing and thawing that should not exceed 10 times. After opening, within 24 hours for short-term preservation, keep sealed at $2-8^{\circ}\text{C}$ in dark. For long-term preservation after opening, keep the lid sealed at $-20^{\circ}\text{C}\pm 5^{\circ}\text{C}$ away from light.

Applicable instruments

Fluorescence microscopy imaging systems, including fluorescence microscopy and filter sets suitable for DAPI (367/452), Green (495/517), and Orange (547/565).

Sample requirements

1. Applicable specimen type: paraffin embedded specimen of surgical resection or biopsy tissue.
2. The tissue in vitro should be fixed with 4% neutral formaldehyde fixative within 1 hour. After the tissue is fixed, it is often dehydrated and embedded in paraffin.

Pretreatment

Baking: Slides heating at 80°C for 30min or 65°C for 2h or overnight.

Dewaxing: According to the customer laboratory protocol (Commonly with Xylene for 15min).

Hydration: Take out the slides and put them respectively into 100%, 85% and 70% EtOH at room temperature for 3 minutes each.

Your authentic partner in molecular cytogenetics

Take out the slides, and immerse them in deionized water for 3 minutes. Remove the excess of water on the slides by air-drying. Permeation: Immerse the slides in deionized water at 100°C and boil continuously for 20-40 minutes (Conventional 20min). Remove the excess of water on the slides by air-drying.

Digestion: Protease enzymic digestion at 37°C for 10-40 minutes. Mix the protease work buffer (50mmol HCl) and the 10x protease solution (Pepsin concentration 0.5%) in a proportion of 9:1 to prepare the enzymatic digestion solution.

Washing: Wash with 2xSSC at room temperature for 5 minutes.

Dehydration: Take out the slides and dehydrate in 70%, 85%, and 100% gradient ethanol at room temperature for 2 minutes each time. Remove the excess of EtOH solution on the slides by air-drying.

Denaturation or hybridization

The following operations should be performed in a darkroom.

1. Take out the probe, let it stand at room temperature for 5min, turn it upside down with force, mix the probe well, centrifuge it briefly (do not vibrate with vortex apparatus), drop 10µl into the hybridization area of the cell drop, cover the 22mm×22mm cover glass immediately, the probe should be evenly spread under the cover glass without bubbles, and seal the edge with rubber (the edge sealing must be thorough to prevent the dry slide from affecting the test results in the hybridization process).
2. Place the glass slide in the hybridization instrument, denature at 85°C for 5 min (the hybridizer should be preheated to 85°C) and hybridized at 42°C for 2-16h.

Washing

The following operations should be performed in a darkroom.

1. Take out the hybridized glass slides, remove the rubber on the coverslip and immediately immerse the slides in a 2xSSC solution for 5 seconds and remove the coverslip.
2. Place the slides in a 2×SSC at room temperature for 1 min.
3. Take out the slides and immerse in a preheated at 68°C 0.3% NP-40/0.4xSSC (Preparation of 0.3% NP-40/0.4xSSC: For 1L preparation, take 3mL NP-40 and 20mL 20xSSC, dissolve fully, mix well, and use 1M NaOH to adjust the pH to 7.2). solution and wash for 2min.

- Remove the slides and immerse in a 37°C preheated deionized water, wash for 1 min and dry the slides naturally in the dark.

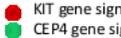
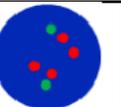
Counterstaining

The following operations should be performed in a darkroom.

10µL DAPI compound dye is dropped in the hybridization area of the glass slide and immediately covered. The suitable filter is selected for glass slide observation under the fluorescence microscope.

FISH results observation

Place the stained slides under a fluorescence microscope and confirm the cells area under a low magnification objective (10×). Under magnification objective (40×) a uniform cells distribution is observed. Then the nuclei FISH results are observed under the high magnification objective (100x).

	
	Negative: 2 orange 2 green
	Positive : n orange 2 green(n≥3)

Precautions

- Please read this manual carefully before testing. The testing personnel shall receive professional technical training. The signal counting personnel must be able to observe and distinguish orange red and green signals.
- When testing clinical samples, if it is difficult to count the hybridization signals and the samples are not enough to repeat the retest, the test will not provide any test results. If the amount of cells is insufficient for analysis, again, the test will not provide test results.
- The formamide and DAPI counterstaining agent used in this experiment have potential toxicity or carcinogenicity, so they need to be operated in the fume hood and wear masks and gloves to avoid direct contact.
- The results of this kit will be affected by various factors of the sample itself, but also limited by enzyme digestion time, hybridization temperature and time, operating environment and limitations of current molecular biology technology, which may lead

Your authentic partner in molecular cytogenetics

to wrong results. The user must understand the potential errors and accuracy limitations that may exist in the detection process.

5. All chemicals are potentially dangerous. Avoid direct contact. Used kits are clinical wastes and should be properly disposed of.
6. This product is for clinical diagnosis and scientific research.



[Manuscript version and approval date]

Manual version: [V1.0 revised on 21 September 2022](#)

Approval date: [15 March 2022](#)